



Handling



Dispensing



Sieving



Blending




Container



Cleaning

## Blending in Pure Culture

Safe, simple, compact blending solutions

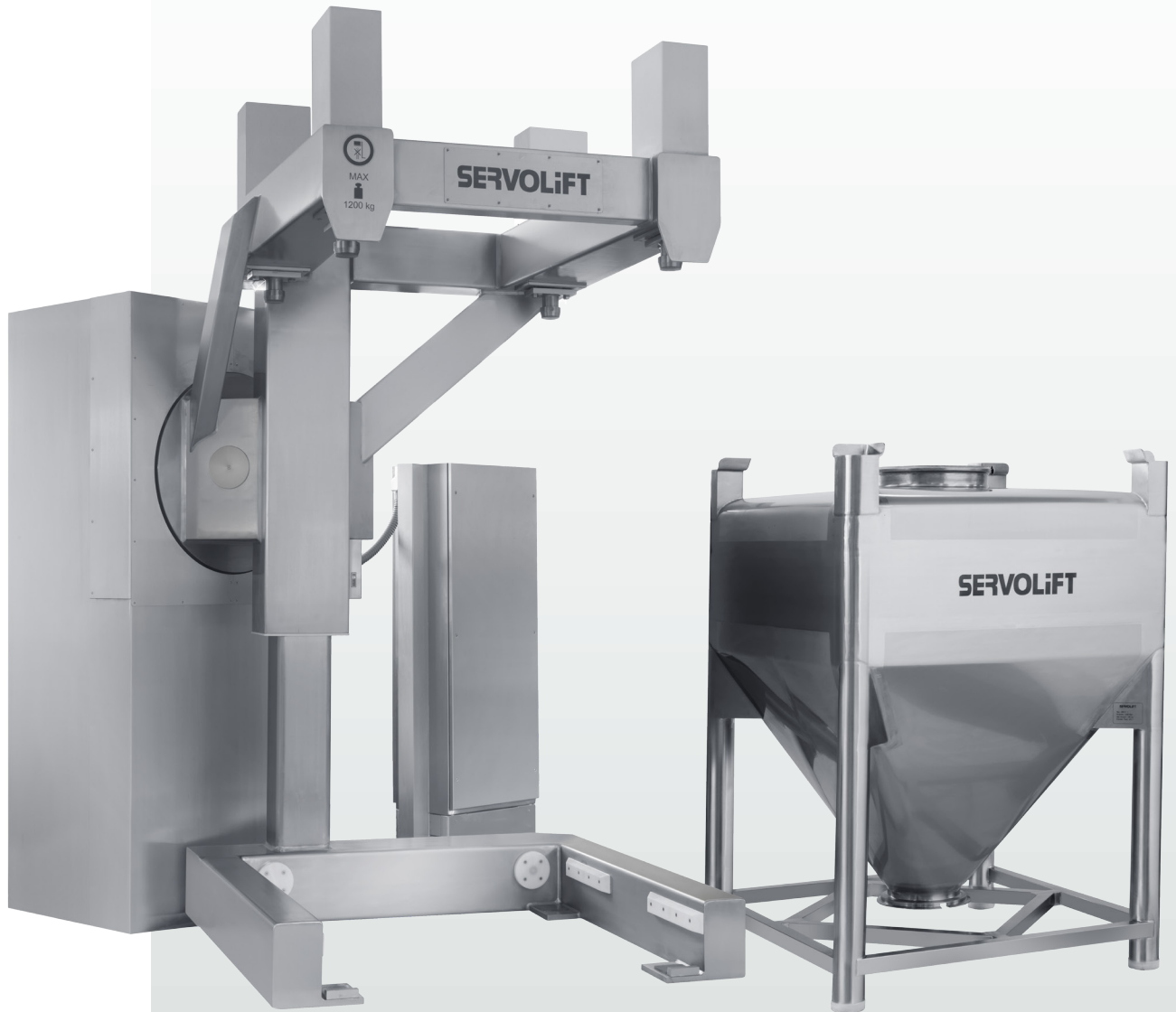

**SERVOLIFT**  
*lifetime solutions*

Your specialist for  
complete handling-, blending-  
and cleaning technologies

SERVOLIFT GmbH  
Albert-Einstein-Straße 9  
77656 Offenburg Germany

T. +49 (0) 781 6100 0  
F. +49 (0) 781 6100 400

info@servolift.de  
www.servolift.de

### ► Blending with a vertical clamping system

This version is a useful alternative to the Servolift standard container blender, especially for containers without clamp bars. The supporting arm is designed so that the container can be positioned by means of a forklift.

The container is positioned on the lower platform, where a jack bar presses it against four compensation cylinders mounted in the upper rigid part of the supporting arm. Thanks to these cylinders, any container tolerances are compensated for and the containers are securely clamped.



Container	Container with stacking corners
Max Load capacity	Up to 2500 kg, according to customer specifications
Safety	Safety light curtain, laser scanner
Installation	Stand-alone, in-wall installation
The following accessories can be installed:	Spraying device for adding liquids, Chopper, NIR spectroscopy
Additional features	can be combined with a CIP system FDR CFR 21 Part 11, FDA-approved construction available
Qualification	IQ, OQ, DQ
Control systems	Siemens S7, Allen-Bradley